


## Matbox®II

### A silicone carrier to protect your surgical instruments

AREX® has designed a one-piece silicone mat sterilization container, to carry and keep safe your surgical instruments.

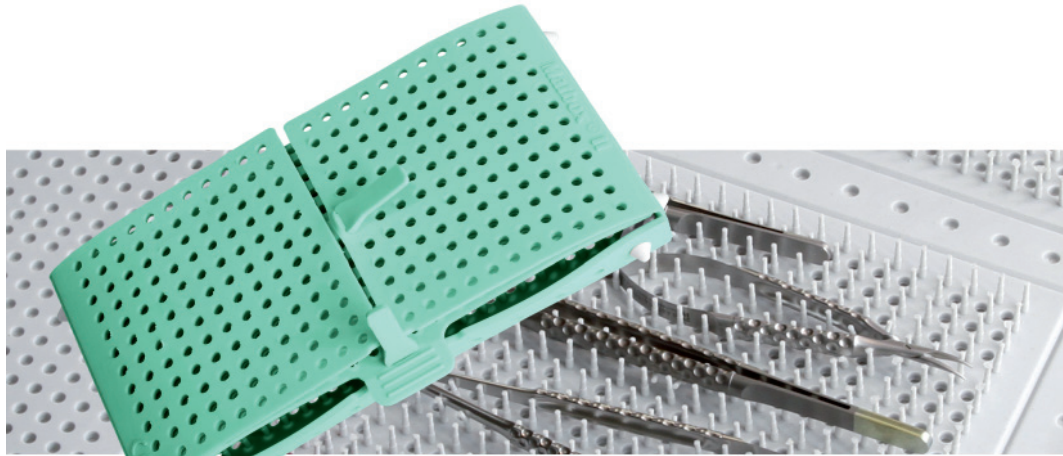

#### TECHNICAL DATA

Platinum cured silicone elastomer  
Silicone hardness: 70 Shore A  
Silicone elastomer complying with USP VI & ISO 10993  
Colouring agent: less than 1% (FDA approved)  
Usable dimensions: 100 x 195 x 20 mm  
Overall dimensions: 115 x 215 x 28 mm

#### AVAILABLE IN A WIDE RANGE OF COLOURS



ref.: BN A OR RO RF GR VT VA NO BL

PROD. REF.: MBX, FOLLOWED BY TWO LETTERS SPECIFYING COLOUR CODE.



MANUFACTURER:  
STERNE SAS  
ZAC du MIN - rue Jean Monnet  
BP50042  
F84302 CAVAILLON Cedex2 - France

[www.arex.fr](http://www.arex.fr)

Design-photo: e-monceaux.fr / Réf. 224 éd.01 - © AREX® - 2013


DISTRIBUTED BY

## Matbox®II

### Effects of silicone spikes design

#### MATBOX®II BASE

Rows of pliable silicone spikes, opposite lateral openings, preventing instrument tips from protruding out of the Matbox®II.

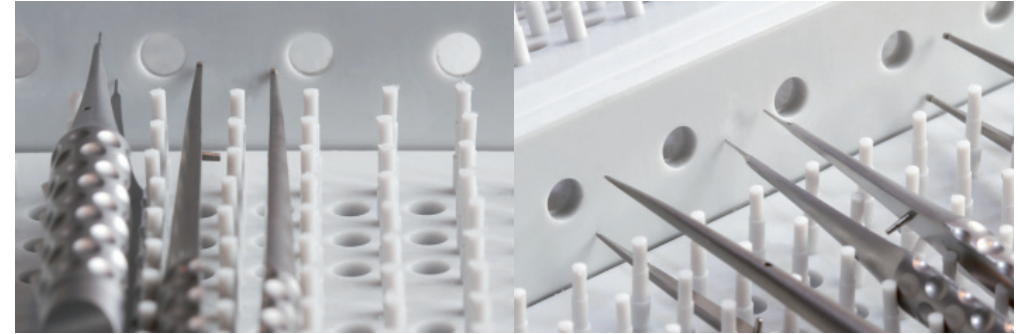

#### MATBOX®II TOP

Rows of pliable silicone spikes staggered with spikes in Matbox®II base, thus gently restraining instruments.

